

### VIS AND UV/VIS SPECTROPHOTOMETERS



# New generation of VIS and UV/Vis spectrophotometers color touch screen display, spectrum scan on-board<sup>1</sup>, USB flash drive recovery

#### PASSION - EXPERIENCE - PERFORMANCE - COMFORT - CONVENIENCE



#### Smart design

A modern design, simple, compact, with all the most common features available directly on the interface. ONDA Touch Series is the perfect instrument for your lab bench top in term of space and solidity




#### **Touchscreen**

A 5" TFT high definition color display for a new level of experience Easy interface, build-in features (photometry, quantitation, spectrum scan<sup>1</sup>) In addition to File management, System Utility and Performance verification



#### **Connectivity and Storage**

USB Flash drive connection for data and methods recovery PC remote control<sup>2</sup>, RS-232 port for micro printer connection, Internal memory to save your methods and data

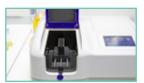

#### Easy access, easy maintenance

Removable front panel, anticorrosive material holder, leaky hole inside sample compartment: all designed to protect the instrument from accidental spills, time wear and damage



#### Lamp room

Removing the top panel for easier and faster access to the source compartment, allows directly control, check of alignment and lamp replacement



#### Multi language

Select the language of your preference and improve your experience with Onda Touch instruments

English, German, Spanish, French, Italian, Portuguese

#### **Applications**

Winery



**Brewery** 



Process water



Waste water



Vegetable Oil/fat

Drinking water



Food&Beverage



**Bottled** water



Sugar industry



Bacterial cell number

<sup>&</sup>lt;sup>1</sup> Spectrum feature available only on SCAN models

<sup>&</sup>lt;sup>2</sup> optional software required

#### THE NEXT GENERATION OF VIS AND UV/VIS SPECTROPHOTOMETERS

#### **New User Interface**

With the new user interface the ONDA Touch series spectrophotometers are ideal for laboratory routine quality control, excellent for applications in education. Work with these innovative instruments as with your smartphone or tablet and enjoy the intuitive, icon-based menu.





#### **Photometry**

Basic photometric measure: Absorbance Transmittance Energy



#### File Management

Save, delete or copy data and methods on internal memory and USB flash drive. Import/export to CSV and TXT files



#### Quantitation

Creation of calibration curve from standard solutions Comfortable calibration of customized methods



#### System utility

System calibration, lamp management, date/time, memory, multi-language



#### Spectrum Scans<sup>1</sup>

High-precision live scans Analysis and automatic peak detection



#### **Performance verification**

Regularly check the instrument functionality for the quality assurance of your results



Internal memory



USB flash drive



Data transfer USB ↔ device



Export to CSV file



Export to TXT file



Data point/Peak detection



Print data and methods



Calibration Curve




Method Detail

#### **Calibration Report**

All Onda Spectrophotometers are supplied with a *Calibration Report* executed according to ISO 17025 requirements, for measuring with traceability.

On the **Calibration Report** are reported:

- Wavelength accuracy
- Photometric accuracy
- Stray light test
- Measurement uncertainty
- Traceability statement

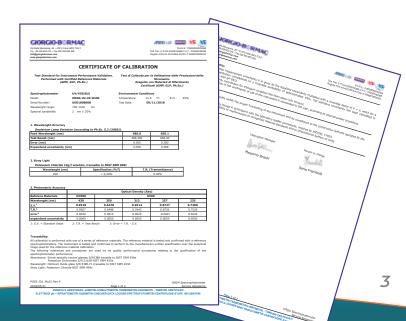

















### **UV-31 SCAN**

**UV/Visible Spectrophotometer** 190 - 1100 nm 2 nm bandwidth Scan speed 4200nm/min















UV-31 Scan is the complete scanning UV-Vis spectrophotometer, ideal for high level educational, the best choice in precision for your qualitative and quantitative analysis. High level of accuracy for sensible detections. The free software included *Easy Uv Basic* allows to perform quantitative, kinetic and photometric single wavelength analysis.

- Low level of stray light for high linearity performance
- Photometry, quantitation, spectrum scan
- Self-calibration system at each switching on
- Preheating countdown
- USB flash drive for data and methods recovery





























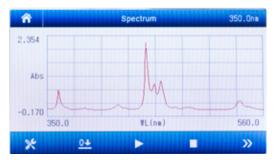

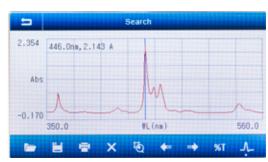

Spectrum scan and detection of peaks

Part No.	Description
11000082	ONDA TOUCH UV-31 SCAN Spectrophotometer, TOUCH SCREEN color display, range 190 - 1100 nm, bandwidth 2 nm, single beam, flash memory USB port, multilanguage OS. Functions: photometric, quantitation, SPECTRUM scan, file management. Including 4 optical glass cells, 2 quarz cells, 10mm 4-cell holder, power cable, dust cover, USB flash drive and Calibration Report, software Easy UV Basic.











### V-11 SCAN

Visible Spectrophotometer 320 - 1100 nm 4 nm bandwidth Scan speed 4200nm/min















*V-11 Scan* is the scanning Visible spectrophotometer, ideal for high level colorimetric analysis, quantitative detections and qualitative control in the visible range.

The free software included *Easy Uv Basic* allows to perform quantitative, kinetic and photometric single wavelength analysis.

- Photometry, quantitation, spectrum scan
- Self-calibration system at each switching on
- Preheating countdown
- USB flash drive for data and methods recovery

































Easy calibration with the online instructions

Part No.	Description
11000062	ONDA TOUCH V-11 SCAN Spectrophotometer, TOUCH SCREEN color display, range 320 - 1100 nm, bandwidth 4 nm, single beam, flash memory USB port, multilanguage OS. Functions: photometric, quantitation, SPECTRUM scan, file management. Including 4 optical glass cells, 10mm 4-cell holder, power cable, dust cover, USB flash drive and Calibration Report, software Easy UV Basic.

















**UV-21** 

**UV/Visible Spectrophotometer** 195 – 1050 nm 4 nm bandwidth













UV-21 is the UV-Vis spectrophotometer, ideal for routine analysis and quantitative detections. The functionality of basic photometry and quantitation make it the perfect instruments for fast and accurate analysis in many application fields. The free software included *Easy Uv Basic* allows to perform quantitative, kinetic and photometric single wavelength analysis.



























P	aı	rt	٨	lo.

#### Description

11000072

ONDA TOUCH UV-21 Spectrophotometer, TOUCH SCREEN color display, range 195 - 1050 nm, bandwidth 4 nm, single beam, flash memory USB port, multilanguage OS. Functions: photometric, quantitation, file management. Including 4 optical glass cells, 2 quarz cells, 10mm 4-cell holder, power cable, dust cover, USB flash drive and Calibration Report, software Easy UV Basic.









Visible Spectrophotometer 325 - 1000 nm 4 nm bandwidth



















V-10 Plus is the Visible spectrophotometer for colorimetric application and routine measurements, widely used in basic photometry, **OD600** in microbiological application and wine color analysis. Solid design and high quality components ensures excellent performance.

ONDA V-10 Plus is supplied with Basic Software for PC (Windows based OS) for the remote control of the device; the application range is expanded to kinetic measurement, improving data analysis and file management.

In combination with the optional ONDA ProWine Color Software for PC it is the best aid for fast and accurate color analysis of wine (Chromatic Characteristics OIV-MA-AS2-07B)

Part No.	Description
11000052	ONDA V-10 Plus Spectrophotometer, including: 4 optical glass cells, 4-position cell holder, power cable with EU
11000052	plug, Basic PC software, dust cover, <i>Calibration Report</i>

### **Cuvettes**

Optical glass 320 – 1400 nm						
Path length (mm)	1	10	20	30	50	100
External dimension WxLxH (mm)	12.5x3.5x45	12.5x12.5x45	12.5x22.5x45	12.5x32.5x45	12.5x52.5x45	12.5x102.5x45
Volume (ml)	0.4	3.5	7.0	10.5	17.5	35.0
Polished windows	2	2	2	2	2	2
Teflon lid	<b>√</b>	<b>√</b>	<b>√</b>	<b>√</b>	<b>√</b>	✓
Number of cells	2	4	4	4	4	4
Part number	11000382	11000252	11000262	11000272	11000282	11000292

Quartz 190 – 1400 nm						
Path length (mm)	1	10	20	30	50	100
External dimension WxLxH (mm)	12.5x3.5x45	12.5x12.5x45	12.5x22.5x45	12.5x32.5x45	12.5x52.5x45	12.5x102.5x45
Volume (ml)	0.4	3.5	7.0	10.5	17.5	35.0
Polished windows	2	2	2	2	2	2
Teflon lid	<b>√</b>	<b>√</b>	<b>√</b>	<b>√</b>	<b>√</b>	<b>√</b>
Number of cells	2	2	2	2	2	2
Part number	11000522	11000302	11000312	11000322	11000332	11000342

Special cells						
Path length (mm)	10	10	10	10	10	10
External dimension WxLxH (mm)	12.5x12.5x45	12.5x12.5x45	12.5x12.5x45	12.5x12.5x45	12.5x12.5x45	12.5x12.5x45
Beam height (mm)	15	15	15	15	15	15
Volume (μl)	100	200	500	1000 (1ml)	700	700
Material	Quartz	Quartz	Quartz	Quartz	Optical glass	Quartz
Range (nm)	190 -1400	190 -1400	190 -1400	190 -1400	320 -1400	190 -1400
Flowing flux	_	_	_	_	<b>√</b>	<b>√</b>
Polished windows	2	2	2	2	2	2
Teflon lid	<b>√</b>	<b>√</b>	<b>√</b>	<b>√</b>	_	_
Number of cells	1	1	1	1	1	1
Part number	11000352	11000362	11000372	11000392	11000492	11000502

## In the following pages useful guide for choice of holder and accessories suitable for your spectrophotometer is provided.

#### **Software Pro-Wine Color for V-10 Plus**

**Pro-Wine Color** is the application software based on Windows operating system for the remote control of **ONDA V-10 Plus Spectrophotometer.**The special application measures absorbance at 420nm, 520nm and 620nm, and calculates automatically Color Intensity and Shade, considering also a Dilution Factor.

# 141

#### Functions:

Photometry, Quantitative, *Kinetics, Multi wavelength, Wine Chromatic Analysis* (Color Intensity and Shade compliance to OIV-MA-AS2-07B)

#### Code to order:

**11000962**, including USB cable, USB dongle key, User manual (the connection PC-instrument requires two USB port on PC)

	4 position for up to 10mm path	4 position for up to 50mm path	4 position for up to 100mm path	Single micro cell, 10mm path	Single test tube ø 8-24mm (max height 110mm)	5° Reflectance measurement	Solid sample (50x30x2mm)	8 position, automatic, 10mm path	5 position automatic Up to 100mm path	Test tube ø 16 mm
						30	_			4
UV-31 SCAN										
3	<b>/</b>	<b>/</b>	<b>/</b>	<b>/</b>	-	<b>\</b>	<b>/</b>	<b>/</b>	<b>\</b>	<b>/</b>
UV-21										
(F)	<b>\</b>	<b>/</b>	<b>\</b>	<b>/</b>	-	<b>\</b>	<b>/</b>	<b>\</b>	<b>/</b>	<b>\</b>
V-11 SCAN										
3	<b>/</b>	<b>/</b>	<b>/</b>	<b>/</b>	-	<b>\</b>	<b>/</b>	<b>/</b>	<b>/</b>	<b>/</b>
V-10 PLUS										
	<b>\</b>	<b>\</b>	<b>\</b>	<b>\</b>	<b>√</b>	<b>√</b>	<b>\</b>	-	-	-
Part number	11000202	11000212	11000222	11000232	11000402	11000412	11000422	11000242	11000432	11000532



#### IQ/OQ paper manual

Printed IQOQ manual, to verify the installation and operation of the spectrophotometer on site.

Part number: 11000192

#### **Software Easy UV for ONDA TOUCH Series**

**Easy UV** is the application software based on Windows operating system for the remote control of **ONDA Touch Series instruments.** The application range, data analysis and file management function are greatly expanded, for example to clinical application (hemoglobin oxygen saturation) and basic biotechnology research (range DNA 0.2-70  $\mu$ g/ml).



#### Functions:

Photometry, Quantitative, Kinetics, *Multi wavelength, Spectrum scan*, *DNA/protein* (UV methods)







#### Code to order:

11000972, including USB Installer, USB cable, USB dongle key, User manual (the connection PC-instrument requires two USB port on PC)

A-1101 Sipper (flowing system)	A-1201 Peltier (Thermostatic system)	A-100 1Sipper/ Peltier (Thermostat flowing system)	Micro Printer (normal paper)	Thermal Printer	PC Software Pro-Wine Color	PC Software ONDA Touch Easy UV	Tungsten- Halogen Iamp (VIS) Vertical filament	Tungsten- Halogen Iamp (VIS) Horizontal filament	Deuterium lamp (UV)
$\checkmark$	<b>/</b>	<b>/</b>	<b>✓</b>	-	-	<b>✓</b>	-	<b>✓</b>	<b>✓</b>
$\checkmark$	<b>/</b>	<b>/</b>	<b>/</b>	-	-	<b>✓</b>	-	<b>✓</b>	<b>✓</b>
<b>/</b>	<b>✓</b>	<b>✓</b>	<b>/</b>	-	-	<b>/</b>	-	<b>✓</b>	-
<b>/</b>	-	-	-	<b>\</b>	<b>✓</b>	-	<b>/</b>	-	-
11000552	11000562	11000572	11000442	11000452	11000962	11000972	11000912	11000942	11000922

## In the following pages useful guide for choice of the appropriate holder suitable for different applications.

Cuvette 10x10 n cod.1100	4 position for up to 10mm path (standard)	4 position for up to 50mm path	4 position for up to 100mm path	Single micro cell, 10mm path
1mm path, 12.5x3.5x45mm	cod.11000582 required	-	-	-
10mm path, 12.5x12.5x45mm	<b>✓</b>	<b>✓</b>	<b>✓</b>	-
20mm path, 12.5x22.5x45mm	-	<b>✓</b>	$\checkmark$	-
30mm path, 12.5x32.5x45mm	-	$\checkmark$	$\checkmark$	-
50mm path, 12.5x52.5x45mm	-	<b>✓</b>	$\checkmark$	-
100mm path, 12.5x102.5x45mm	-	-	<b>✓</b>	-
Test tube ø16mm, 90mm height	-	-	-	-
Solid sample 50x30x2mm	-	-	-	-
Solid sample 20x20x4mm	-	-	-	-
Micro-cell 12.5x12.5x45mm	-	-	-	$\checkmark$
Flow-cell 12.5x12.5x45mm	-	-	-	$\checkmark$
Part Number	11000202	11000212	11000222	11000232

Single test tube ø 8-24mm (max height 110mm)	Test tube ø 16 mm	5° Reflectance attachement	Solid sample (50x30x2mm	8 position, auto- matic, 10mm path	5 position, automatic, up to 100mm path
*	•		1		
-	-	-	-	cod.11000582 required	cod.11000582 required
-	-	-	-	<b>✓</b>	<b>✓</b>
-	-	-	-	-	<b>✓</b>
-	-	-	-	-	<b>✓</b>
-	-	-	-	-	<b>/</b>
-	-	-	-	-	<b>/</b>
<b>/</b>	<b>/</b>	-	-	-	-
-	-	-	<b>/</b>	-	-
-	-	<b>/</b>	-	-	-
-	-	-	-	-	-
-	-	-	-	-	-
11000402	11000532	11000412	11000422	11000242	11000432

Technical Specification	UV-31 SCAN	UV-21	V-11 SCAN	V-10 PLUS		
Spectrophotometer		UV-VIS		VIS		
Optical system		Sir	ngle beam			
Light source	Tungsten-Halog	Tungsten-Halogen / Deuterium Lamps Tungs				
Detector						
Spectral bandwidth (nm)	2	4	4	4		
Wavelenght range (nm)	190 - 1100	195 – 1050	320 – 1100	325 - 1000		
Wavelenght accuracy (nm)	±0.8	±1.0	±1.0	±2		
Wavelenght repeatability (nm)	±0.2	±0.4	±0.2	±0.8		
Wavelenght resolution (nm)	0.1	0.1	0.1	0.1		
Photometric range		-0.3 to 3 Abs / 0 to 2	200 %T / 0 to 9999.9 C	Conc		
Photometric accuracy		±0.5 %T	at 0 - 100 %T			
Photometric repeatability		±0.2 %T at 0 - 100 %	6T	±0.3 %T at 0 – 100 %T		
Photometric resolution (Abs)			0.001			
Stability		±0.002 Abs/hr (at 50	00nm after 2hours warn	nup)		
Stray light	<0.1 %T	<0.2 %T	<0.1 %T	<0.3 %T		
Display	5 inches co	olor touchscreen display, f	ixed, 480x272 px	LCD display 128x64 px		
Keypad		touchscreen		Membrane keypad		
Standard cell holder		4 positions,10mm pat	h square cells, manual l	holder		
Sample compartment		Accessible from top and	front	Accessible from top		
Basic photometry	✓	✓	✓	✓		
Quantitation	✓	✓	✓	✓		
Spectrum scan	✓	-	✓	-		
File management	✓	✓	✓	-		
System management	✓	✓	✓	-		
Performance verification	✓	✓	✓	-		
Connectivity		USB-A (flash drive) USB-B (PC connectio RS-232 (printer) <sup>2</sup>	n)¹	USB-B (PC connection) Serial port (printer) <sup>2</sup>		
Memory	E	236 kB (internal memo expandable with USB flast		200 raw of data		
Languages	English, Sp	anish, French, German, It	alian, Portuguese	English		
Dimensions WxDxH (mm)	456 x 360 x 185		490 x 360 x 210			
Weigth (kg)	10.7 10.5			12		
Power supply		100 – 240	V AC, 50/60 Hz			
Power consumption (W)		120		75		
Warranty			2 years <sup>4</sup>			
Part number	11000082	11000072	11000062	11000052		

<sup>&</sup>lt;sup>1</sup> Optional PC software required

DISTRIBUTED BY:

GUARANTEED BY:

GIORGIO-BORMAC





<sup>&</sup>lt;sup>2</sup> Optional thermal printer required

<sup>&</sup>lt;sup>3</sup> USB flash drive included with ONDA Touch series

 $<sup>^{\</sup>mbox{\tiny 4}}$  On electronic and optical part (light source excluded)